

SIEMENS

FLUOROSPOT H/HK

AX

Maintenance Protocol

Customer:

Department:

Room:

Address:

Contact person:

Telephone:

Cust. specific no.:

Customer no:

The Maintenance Instructions
RX41-020.101.04.01.02
are required for this protocol

© Siemens AG 1997

The reproduction, transmission or use of this document or its contents is not permitted without express written authority. Offenders will be liable for damages. All rights, including rights created by patent grant or registration of a utility model or design, are reserved.

Register 9

Print Number: RX41-020.105.04.01.02

Replaces: n. a.

English

Doc. Gen.Date: 05.97

SIEMENS office:

Address:

Region:

Country:

Contact person:

Telephone:

FSE in charge:

Telephone:

Information regarding the maintenance protocol

The maintenance protocol serves as a quality record of the maintenance performed

Maintenance must be performed in the prescribed intervals.

The maintenance results must be entered in the maintenance protocol.

The page numbers in the check list at the end of the protocol refer to the corresponding pages in the associated maintenance instructions (refer to the cover sheet).

The maintenance protocol must be completely filled out by the FSE in charge.

In other words,

- All fields must be completed. If a field does not apply to the system or if there is no information to be entered, enter 'n.a.' in the field.
- The customer number must be entered in the header of every page (CS No.) so that every page is assigned to a customer.
- In case of complaints, enter the product identification (IVK = WE) for the component as well as the type of complaint in the 'Open items' table provided. Record the open items in the table with the date and signature. Correction of these open points should also be documented in this table along with the date and signature. If there are no open points, draw a line through the entire table and enter the date/initials.
- Enter the values measured during the maintenance in the space / table provided.
- After completing the maintenance, fill out page 3 of the protocol and sign it.

Processing and archiving the maintenance protocol:

The maintenance protocol is considered a record and must be archived. It must be filed after completion of maintenance in the corresponding Register of the Maintenance binder. If necessary, a copy should be distributed to the customer.

System Status:

System:

Ser. No.:

Software Version:

Maintenance contract no.:

Type of contract:

| | |
|---|--------------------------|
| The equipment has no problems Results of the image quality test showed no deviations from the required reference values. | <input type="checkbox"/> |
| The equipment has minor problems that do not restrict its use. However they should be corrected as a preventive measure. Results of the image quality test showed no deviations from the required reference values. | <input type="checkbox"/> |
| The equipment has major problems. For safety reasons, the equipment may not be used until the problem has been resolved. | <input type="checkbox"/> |

| | |
|-----------------|--|
| Location | |
| Date | |
| Name of the FSE | |
| Signature | |

Performed tasks

[illegible]

Open points

| IVK (WE) | Component | Open points | Solved (Date) | FSE (Date) |
|-------------|-----------|-------------|------------------|---------------|
| | | | | |
| | | | | |
| | | | | |
| | | | | |
| | | | | |
| | | | | |
| | | | | |
| | | | | |
| | | | | |
| | | | | |

Explanation of the acronyms:

| Abbrev. | Explanation |
|---------|---|
| SI | Safety Inspection |
| SIE | Safety Inspection Electrical Safety |
| SIM | Safety Inspection Mechanical Safety |
| PM | Preventive Maintenance |
| PMP | Preventive Maintenance Preventive Parts Exchange, External Inspection, etc. |
| PMA | Preventive Maintenance Adjustments |
| PMF | Preventive Maintenance Function, Operating-Value Check |
| Q | Quality Check |
| QIQ | Quality Check Image |
| QSQ | Quality Check System |
| SW | Software Maintenance |

Measurement devices

Please enter all measuring devices used for maintenance in the table.

| Measurement Device | Type | Serial-No. | Date of use | Calibrated until |
|--------------------------|------|------------|-------------|------------------|
| TV Dynamics Test | | | | |
| Luminous intensity meter | | | | |
| | | | | |
| | | | | |
| | | | | |
| | | | | |
| | | | | |
| | | | | |
| | | | | |
| | | | | |
| | | | | |
| | | | | |
| | | | | |
| | | | | |
| | | | | |
| | | | | |
| | | | | |

Image Quality Test

Indirect Dose Rate Control (IQ Chap. 4.1.2)

Full format, **all** control steps

Prefiltering: 2.1 mm precision radiation filter

| Operating Conditions | | | | Initial Values | | | | Mainten. |
|-------------------------------------|--------------------------|-------|------|--------------------|------|--------------------|--|----------|
| Control step | Fluoro Progr. (nGy/s) | Curve | Zoom | + | | - | | |
| | | | | - 1 kV - 10% mA | Act. | + 1 kV + 10% mA | | |
| Step 1 Fluoro 1 Step 2 | | | | kV mA | | | | |
| | | | | kV mA | | | | |
| Step 1 Fluoro 2 Step 2 | | | | kV mA | | | | |
| | | | | kV mA | | | | |
| Step 1 Fluoro 3 Step 2 | | | | kV mA | | | | |
| | | | | kV mA | | | | |

Zoom formats, Fluoro 1, Step 1

| I.I. Format / I.I. | Initial Values | | | | Mainten. |
|--------------------|--------------------|-----|--------------------|--|----------|
| | + | | - | | |
| | - 1 kV - 10% mA | Act | + 1 kV + 10% mA | | |
| Zoom 1 | kV mA | | | | |
| Zoom 2 | kV mA | | | | |
| Zoom 3 | kV mA | | | | |

Remarks: _____

Indirect Dose Control (IQ Chap. 8.1.) DR

| Dose | kV | F/s | Scene | Focus | I.I. For- mat (cm) | Zoom Dose fac. | Initial Values Qg/n | | | Mainten. |
|------|----|-----|-------|-------------------------------------|-----------------------|-------------------|---------------------|-----|-------|----------|
| | | | | | | | - | Act | + | |
| | | | | | | | - 15% | | + 15% | |
| 100 | 70 | 1 | 10s | <input checked="" type="checkbox"/> | | | | | | |
| | | | | | | | | | | |
| | | | | | | | | | | |
| | | | | | | | | | | |

DSA

| Dose | kV | F/s | Scene | Focus | I.I. For- mat (cm) | Zoom Dose Fac. | Initial Values Qg/n | | | Mainten. |
|------|----|-----|-------|-------------------------------------|-----------------------|-------------------|---------------------|-----|-------|----------|
| | | | | | | | - | Act | + | |
| | | | | | | | - 15% | | + 15% | |
| 500 | 70 | 0.5 | 10s | <input checked="" type="checkbox"/> | | | | | | |

Remarks: _____

Monitor Adjustment (IQ Chap.8.3)

DFR Test Images

| Monitor | | | Initial Values (cd/m ²) | | | Mainten. |
|--------------------------|------|-------|--|---|-----|---|
| Luminous intensity | | | - | Act | + | |
| Examination room | Live | black | 0.7 | | 0.9 | |
| | | white | 240 | | 280 | |
| | Ref | black | 0.7 | | 0.9 | |
| | | white | 240 | | 280 | |
| 5% / 95% fields visible? | | | | <input type="checkbox"/> yes <input type="checkbox"/> no | | <input type="checkbox"/> yes <input type="checkbox"/> no |
| Control room | Live | black | 0.7 | | 0.9 | |
| | | white | 240 | | 280 | |
| | Ref | black | 0.7 | | 0.9 | |
| | | white | 240 | | 280 | |
| 5% / 95% fields visible? | | | | <input type="checkbox"/> yes <input type="checkbox"/> no | | <input type="checkbox"/> yes <input type="checkbox"/> no |

NOTICE

On SIMOMED monitors, switch off the ambient light sensor for the measurement!

Remarks: _____

Dynamics Test with Fluoroscopy (IQ Chap. 5.8.1) with VIDEOMED SX**B/BA-Signal Amplitudes**

Focus-I.I. distance (SID) = _____ cm


I.I. format ☐ 22 cm; ☐ 23 cm; ☐ 27 cm; ☐ 28 cm

- with grid
- Dynamics test (close to I.I.)
- without heart contour collimator
- without capillary test
- B-Signal **without bias light**

Fluoro curve _____ \odot = _____ $\mu\text{Gy/s}$

| B-/BA-Signal (mV) | Field No. | Initial Values (mV) | | | Mainten. |
|--------------------------------------|-----------|---------------------|-----|-----|----------|
| | | - | Act | + | |
| B-Signal Input TV CCU | 2L | — | | — | |
| | 1 | 65 | | 95 | |
| | 5R | — | | — | |
| Fluoro dyn. | | | | | |
| BA-Signal Output TV CCU | 2L | ≥ 1300 | | — | |
| | 1 | 300 | | 400 | |
| | 5R | — | | — | |
| Fluoro kV | | — | | — | |

*1 Procedure with SIREGRAPH D2/D3/CF

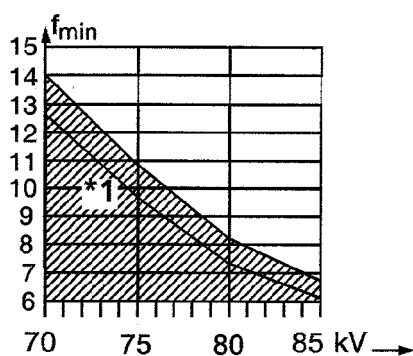

Dynamics factor (f) = level/field 2L: level/field 1

Remarks: _____

Dynamics Test with DR and DSA (IQ Chap.8.5)

CAUTION

Do not select any additional filtering in the collimator

DR: Dose 100, single image

Iris corr. values _____

Focus-I.I. distance (SID) = _____ cm

- I.I. format ☐ 22 cm; ☐ 23 cm; ☐ 27 cm
or ☐ 28 cm

- Dynamics test (close to I.I.)

- with heart contour collimator

- with capillary test

- with grid

upper window value: W1 = 511

lower window value: W2 = 256

- Measure BA-Signal at DFR output

DSA: Dose 200, F/s = 6, scene 5s

Iris corr. values _____

Amplitudes

| | Field No. | Initial Values (mV) | | | Mainten. |
|-----------------------------------|-----------|---------------------|-----|-------|----------|
| | | - | Act | + | |
| DR BA-Signal (mV) | 2L | 520 | | 620 | |
| | 1 | 200 | | 250 | |
| | 5R | — | | — | |
| DSA * BA-Signal (mV) | 2L | — | | < 650 | |
| | 1 | 190 | | 260 | |
| | 5R | — | | — | |
| U _a (kV) | DR | — | 70 | — | |
| | DSA | — | 70 | — | |
| U _{Iris} (V) | DR | — | | — | |
| | DSA | — | | — | |

Reference values applicable for 70 kV

* Evaluate Native image (next to last image);

If 2L field would be > 650 mV, field 1 < 190 mV is admissible.

Remarks: _____

Edge Enhancement (IQ Chap.8.5.2)

DR exposure

| Function configured | Maintenance |
|---------------------|-------------|
| yes | |
| no | |

Remarks: _____


Mean Value Calculation (GGM) (IQ Chap. 8.5.3)

Dig. Fluoro

| Function configured | Maintenance |
|---------------------|-------------|
| yes | |
| no | |

Remarks: _____

Resolution (IQ Chap. 8.7) and Minimum Contrast (IQ Chap. 8.8)

- 17 μm Cu strips
- With grid, 
- Resolution test, Type 41 turned 45° and collimate to the test size
- Insert 1.2 mm Cu in the collimator
- E (%) = 0; temp. filter = 6%
- Fluoro 1 at fluoro (\approx 70 kV)
- Dose step 100 with DR (\approx 70 kV)
- Upper window value: W1 = 511; lower window value: W2 = 256

| | Format | Initial Values (cd/m ²) (LP/mm) | | | Mainten. |
|------------------------|-------------|--|-----|---|----------|
| | | - | Act | + | |
| Min. contrast visible? | | — | yes | | |
| Dig. Fluoro @ 1 | Full format | | | — | |
| | Zoom 1 | | | — | |
| | Zoom 2 | | | — | |
| | Zoom 3 | | | — | |
| Min. contrast visible? | | — | yes | | |
| DR exposure | Full format | | | — | |
| | Zoom 1 | | | — | |
| | Zoom 2 | | | — | |
| | Zoom 3 | | | — | |
| Min. contrast visible? | | — | yes | | |
| Hardcopy | Full format | | | — | |
| | Zoom 1 | | | — | |
| | Zoom 2 | | | — | |
| | Zoom 3 | | | — | |

Remarks: _____

Check of the DSA Device (IQ Chap. 8.9)

Focus-I.I. distance (SID)= _____ cm

I.I. format ☐ 22 cm; ☐ 23 cm; ☐ 27 cm; ☐ 28 cm

Window values for DSA image: brightness W1 = 80; contrast: W2 = 220

- With grid
- E (%) = 0
- Dynamics test with heart contour collimator and capillary test 37 90 180
- Trigger a DSA scene: 6 F/s, Dose 200, scene 5s, 70 kV

| Contrast Sensitivity | | | | Initial Values | | | Maintenance *1 | | |
|----------------------|--------|--|--|----------------|---|----|----------------|---|----|
| | | | | 2L | 1 | 5R | 2L | 1 | 5R |
| 3 mm | 16 % | | | | | | | | |
| | 4 % | | | | | | | | |
| | 0.75 % | | | | | | | | |
| 2 mm | 10 % | | | | | | | | |
| | 3 % | | | | | | | | |
| | 2 % | | | | | | | | |
| 1 mm | 7.5 % | | | | | | | | |
| | 4 % | | | | | | | | |
| | 2.5 % | | | | | | | | |

| Tolerance | | n+2 | n | n - 3 | LOG amplifier *2 | | |
|-----------|------|-----|---|-------|------------------|--|--|
| 3 mm | 16 % | | | | | | |

Subtraction *3

Tol. 0 ± 1

| | | | |
|--|--|--|--|
| | | | |
|--|--|--|--|

*1 mark capillaries **not** visible

*2 Enter gray step of capillary

*3 Enter basic gray value of test field

Remarks: _____

Image Disturbances (Artifacts) (Chap. 8.10)

DR Mode (Chap. 8.10.1)

| Images from MP | Type of disturbance, Artifact | Initial Values | Mainten. *1 |
|----------------|-------------------------------|-----------------|-------------|
| 8.2 | Pixel errors/ RAM | | |
| 8.5.1 | | POST PROCESSING | |
| 8.5.1 | Ghost images | | |
| 8.2 / 8.5.1 | Contouring | | |
| *2 | Horizontal jitter | | |
| 8.2 / 8.5.1 | Island formation | | |
| 8.4.3 | Vignetting | | |

Other disturbances: _____

DSA Mode (IQ Chap. 8.10.2)

| Images from MP | Type of Disturbance, Artifact | Initial Values | Mainten. *1 |
|----------------|--------------------------------------|-----------------|-------------|
| 8.2 | Pixel errors/ RAM | | |
| 8.5.1 | | POST PROCESSING | |
| 8.9 | Artifacts caused by exposure | | |
| 8.9 | Logarithming errors | | |
| 8.9 / *2 | Convergence errors | | |
| 8.9 / *2 | Inhomogeneousness (Hum, Line Errors) | | |
| 8.9 / *2 | Microphony | | |
| 8.9 / 8.2 | Island formation | | |
| 8.9 | Noise patterns | | |

Other disturbances: _____

*1 Enter No. 1, 2, 3


1 = No disturbances, artifacts

2 = Slight disturbances, artifacts

3 = Intolerable disturbances, artifacts

*2 Soldering wire cross on middle of I.I.; 0.6 mm Cu prefilter

Digital Tomography (IQ Chap. 8.11)

- Tomo test object on tabletop
- Tomo height approx. 11.5 cm
(depending on tomo test object)
-  , 0.6 mm Cu
- I.I. format 17 cm (20 cm)
- With grid
- Focus-I.I. distance (SID) = _____ cm
- Exposure voltage = _____ kV; Dose 100
- Tomo profile
- max. tomo angle: _____ °
- min. tomo time: _____ s

| | | Initial Values | Mainten. |
|---------------------------------------|------------|----------------|----------|
| * $\Delta S_{Ha} = S_{Ha} - S_H$ (mm) | | | |
| $\Delta S_{HL} = S_{HL} - S_H$ (mm) | | | |
| Rg (LP/mm) | | | |
| Shape and sweep | / o.k. | | |
| | / not o.k. | | |
| Blur | / o.k. | | |
| | / not o.k. | | |

*) Explanations

S_{Ha} = Displayed tomo height (≤ 5 mm)

S_H = Actual tomo height

S_{HL} = Light line height (≤ 5 mm)

Rg = Measured resolution (≥ 2.0)

Remarks:

This page intentionally left blank.

| | | OK not n.a. OK | Page |
|-------------|---|-------------------|---------------|
| 1 | Requirements | | 1 - 1 |
| 1.1. | General Remarks | | 1 - 1 |
| 1.2. | Required Documents | | 1 - 1 |
| 1.3. | Required Test Equipment and Aids | | 1 - 2 |
| 1.4. | Explanation of Symbols | | 1 - 2 |
| 2 | VIDEOMED SX | | 2 - 1 |
| 2.1. | Safety Inspection | | 2 - 1 |
| SI | Checking the Camera Mount. | ○ ○ ○ | 2 - 1 |
| SI | Checking the Camera Housing | ○ ○ ○ | 2 - 1 |
| SI | Checking the Ground Wire | ○ ○ ○ | 2 - 1 |
| 2.2. | Maintenance | | 2 - 1 |
| PM | Camera Optics | ○ ○ ○ | 2 - 1 |
| PM | CCU | ○ ○ ○ | 2 - 1 |
| 3 | FLUOROSPOT H/HK | | 3 - 1 |
| 3.1. | External Visual Inspection | | 3 - 1 |
| PM | Check the Control Console for Damage. | ○ ○ ○ | 3 - 1 |
| 3.2. | Maintenance | | 3 - 1 |
| PM | Opening the Unit | ○ ○ ○ | 3 - 1 |
| PM | Remove dust that has collected on the inside of the unit | ○ ○ ○ | 3 - 1 |
| 3.3. | Inspection of the internal head vents | | 3 - 1 |
| PM | Check the fan for proper function. | ○ ○ ○ | 3 - 1 |
| 3.4. | Monitors | | 3 - 1 |
| PM | Cleaning the Monitor Screens | ○ ○ ○ | 3 - 1 |
| 3.5. | Function Inspection | | 3 - 2 |
| PM | Checking / Read Out the Error Log | ○ ○ ○ | 3 - 2 |
| 4 | Image Quality Test | | 4 - 1 |
| 4.1. | General Remarks | | 4 - 1 |
| Q | Indirect Dose Rate Control (IQ Chap. 4.1.2) | ○ ○ ○ | 4 - 2 |
| Q | Indirect Dose Control (IQ Chap. 8.1.2) | ○ ○ ○ | 4 - 3 |
| Q | Monitor Adjustment (IQ Chap. 8.3) | ○ ○ ○ | 4 - 5 |
| Q | Dynamics Test (Fluoroscopy) (IQ Chap. 5H, 5.8) | ○ ○ ○ | 4 - 6 |
| Q | Dynamics Test (DFR) (IQ Chap. 8.5). | ○ ○ ○ | 4 - 9 |
| Q | Edge Enhancement (IQ Chap. 8.5.2) | ○ ○ ○ | 4 - 12 |
| Q | Mean Value Evaluation (GGM) (IQ Chap. 8.5.3) | ○ ○ ○ | 4 - 12 |
| Q | Resolution (IQ Chap. 8.7) and Minimum Contrast (IQ Chap. 8.8) | ○ ○ ○ | 4 - 13 |
| Q | Checking the DSA Device (IQ Chap. 8.9) | ○ ○ ○ | 4 - 14 |
| Q | Image Disturbances (Artifacts) (Chap. 8.10) | ○ ○ ○ | 4 - 17 |
| Q | Digital Tomo IQ Chap. 8.11) | ○ ○ ○ | 4 - 20 |
| Q | Check of the Tomo Height | ○ ○ ○ | 4 - 21 |
| Q | Checking the Resolution | ○ ○ ○ | 4 - 23 |
| 4.2. | Final Work Steps | | 4 - 24 |

| | | OK not n.a. OK | Page |
|-------|-----------------------------|-------------------|--------|
| 4.2.1 | Generator Fluoro Curves | | 4 - 24 |
| 4.2.2 | SIMOMED H | | 4 - 24 |
| 4.2.3 | Function Checks | | 4 - 24 |
| 5 | Changes to Previous Version | | 5 - 1 |